




Telephone: (650) 697-3600



## **GOAT ANTI-DESMIN ANTIBODY**

**SKU:** EB09106

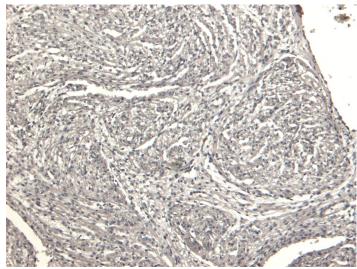

## **SPECIFICATIONS**

Formulation Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

**Unit Size** 100 μg

Storage

Aliquot and store at -20°C. Minimize freezing and thawing. Instructions

Synonym /

intermediate filament protein|FLJ41793|FLJ41013|FLJ39719|FLJ12025|CSM2|CSM1|CMD1I|desmin|DES **Alias** 

**Names** 

**Accession** NP\_001918.3 ID

**Blocking** 

EBP09106 **Peptide** 

Peptide with sequence C-RDGEVVSEATQQQHE, from the C Terminus of the protein sequence according to Immunogen

NP 001918.3.

**Peptide** 

C-RDGEVVSEATQQQHE

Sequence

Purification Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography

Method using the immunizing peptide.

**Shipping** 

Refrigerated Instructions





Email: <a href="mailto:customerservice@vectorlabs.com">customerservice@vectorlabs.com</a>

Telephone: (650) 697-3600

**Predicted** 

Human, Mouse, Dog, Cow

Species Reactive

Species Human, Mouse

Human

Gene ID

Mouse

13346

Gene ID

**Rat Gene ID** 64362

Product Grade

https://prod-vector-labs-pimcore-assets.s3.us-east-1.amazonaws.com/assets/products/image/elite\_medium.png

IHC Results Paraffin embedded Human Heart. Recommended concentration: 4-6μg/ml.

**ELISA** 

Detection

Antibody detection limit dilution 1:128000.

Limit

Approx. 55-60kDa band observed in Human and Mouse Skeletal Muscle lysates and in Human Heart lysates

Western Blot (calculated MW of 53.5kDa according to Human NP\_001918.3, and Mouse NP\_034173.1). Preliminary testing

showed a band at approx 60kDa in Rat and Pig Skeletal muscle lysate. Recommended concentration:

0.01-0.3µg/ml. Primary incubation 1 hour at room temperature.

Application

Type

Pep-ELISA, WB, IHC

## **DOCUMENTS**


• Data Sheet

## **GALLERY IMAGES**



Telephone: (650) 697-3600



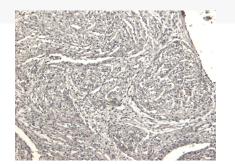







Email: <a href="mailto:customerservice@vectorlabs.com">customerservice@vectorlabs.com</a>

Telephone: (650) 697-3600

