Telephone: (650) 697-3600



 $\textbf{Email:} \ \underline{\textbf{customerservice@vectorlabs.com}}$ 

## GOAT ANTI-TEB4 / MARCH-VI (AA 701 TO 712) ANTIBODY

**SKU:** EB08294



Email: <a href="mailto:customerservice@vectorlabs.com">customerservice@vectorlabs.com</a>

Telephone: <u>(650)</u> 697-3600

250kDa 150kDa 100kDa

75kDa

50kDa

37kDa

25kDa

20kDa

15kDa



Telephone: (650) 697-3600

## **SPECIFICATIONS**

Formulation Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

**Unit Size** 100 μg

Storage

Instructions Aliquot and store at -20°C. Minimize freezing and thawing.

Synonym /

Alias

membrane-associated RING-CH protein VI/RNF176|KIAA0597|membrane-associated ring finger (C3HC4)

6|MARCH6|TEB4|MARCH-VI Names

**Accession** 

ID

NP\_005876.2

**Blocking** 

EBP08294 **Peptide** 

Peptide with sequence C-QGRRVIFQKVKE, from the internal region of the protein sequence according to **Immunogen** 

NP 005876.2.

**Peptide** 

Sequence

C-QGRRVIFQKVKE

Purification Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography

Method using the immunizing peptide.

**Shipping** 

Instructions

Refrigerated

**Predicted** 

**Species** 

Human, Mouse

Reactive

**Species** 

Human

223455

Human

10299 Gene ID

Mouse

**Gene ID** 

**Product** 

Grade

https://prod-vector-labs-pimcore-assets.s3.us-east-1.amazonaws.com/assets/products/image/elite\_medium.png

**ELISA** 

**Detection** 

Antibody detection limit dilution 1:16000.

Limit Western


Approx 100kDa band observed in Human Lung lysates (calculated MW of 103kDa according to NP\_005876.2). Recommended concentration: 0.2-0.5µg/ml. An additional band of unknown identity was also consistently observed at 18kDa. This band was successfully blocked by incubation with the immunising peptide.

**Application** 

**Type** 

Blot

Pep-ELISA, WB

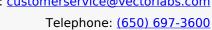



## **DOCUMENTS**

• Data Sheet

## **GALLERY IMAGES**

250kDa 150kDa 100kDa

75kDa

50kDa

37kDa

25kDa

20kDa

15kDa





Email: <a href="mailto:customerservice@vectorlabs.com">customerservice@vectorlabs.com</a>

Telephone: (650) 697-3600